

PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 6 : <b>G01N</b>		A2	(11) International Publication Number: <b>WO 96/13707</b> (43) International Publication Date: <b>9 May 1996 (09.05.96)</b>
(21) International Application Number: <b>PCT/US95/12550</b> (22) International Filing Date: <b>20 October 1995 (20.10.95)</b>		(81) Designated States: CA, JP, European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).  <b>Published</b> <i>Without international search report and to be republished upon receipt of that report.</i>	
(30) Priority Data: 08/326,788 20 October 1994 (20.10.94) US			
(71)(72) Applicant and Inventor: <b>CONNOLLY, James [US/US]; 6131 West 80th Street, Indianapolis, IN 46278 (US).</b>			
(74) Agent: <b>LITZINGER, Jerrold, J.; Suite 600, 4445 Lake Forest Drive, Cincinnati, OH 45242 (US).</b>			
(54) Title: <b>APPARATUS AND METHOD FOR DETERMINING SUBSTANCES CONTAINED IN A BODY FLUID</b>			
<img alt="Diagram of the apparatus showing a test strip 52 being inserted into a device. The device has a display 66 and various buttons and ports labeled 54, 56, 58, 60, 62, 64, 66, 68, 70, 72, 74, 76, 78, 80, 82, 84, 86, 88, 90, 92, 94, 96, 98, 100, 102, 104, 106, 108, 110, 112, 114, 116, 118, 120, 122, 124, 126, 128, 130, 132, 134, 136, 138, 140, 142, 144, 146, 148, 150, 152, 154, 156, 158, 160, 162, 164, 166, 168, 170, 172, 174, 176, 178, 180, 182, 184, 186, 188, 190, 192, 194, 196, 198, 200, 202, 204, 206, 208, 210, 212, 214, 216, 218, 220, 222, 224, 226, 228, 230, 232, 234, 236, 238, 240, 242, 244, 246, 248, 250, 252, 254, 256, 258, 260, 262, 264, 266, 268, 270, 272, 274, 276, 278, 280, 282, 284, 286, 288, 290, 292, 294, 296, 298, 300, 302, 304, 306, 308, 310, 312, 314, 316, 318, 320, 322, 324, 326, 328, 330, 332, 334, 336, 338, 340, 342, 344, 346, 348, 350, 352, 354, 356, 358, 360, 362, 364, 366, 368, 370, 372, 374, 376, 378, 380, 382, 384, 386, 388, 390, 392, 394, 396, 398, 400, 402, 404, 406, 408, 410, 412, 414, 416, 418, 420, 422, 424, 426, 428, 430, 432, 434, 436, 438, 440, 442, 444, 446, 448, 450, 452, 454, 456, 458, 460, 462, 464, 466, 468, 470, 472, 474, 476, 478, 480, 482, 484, 486, 488, 490, 492, 494, 496, 498, 500, 502, 504, 506, 508, 510, 512, 514, 516, 518, 520, 522, 524, 526, 528, 530, 532, 534, 536, 538, 540, 542, 544, 546, 548, 550, 552, 554, 556, 558, 560, 562, 564, 566, 568, 570, 572, 574, 576, 578, 580, 582, 584, 586, 588, 590, 592, 594, 596, 598, 600, 602, 604, 606, 608, 610, 612, 614, 616, 618, 620, 622, 624, 626, 628, 630, 632, 634, 636, 638, 640, 642, 644, 646, 648, 650, 652, 654, 656, 658, 660, 662, 664, 666, 668, 670, 672, 674, 676, 678, 680, 682, 684, 686, 688, 690, 692, 694, 696, 698, 700, 702, 704, 706, 708, 710, 712, 714, 716, 718, 720, 722, 724, 726, 728, 730, 732, 734, 736, 738, 740, 742, 744, 746, 748, 750, 752, 754, 756, 758, 760, 762, 764, 766, 768, 770, 772, 774, 776, 778, 780, 782, 784, 786, 788, 790, 792, 794, 796, 798, 800, 802, 804, 806, 808, 810, 812, 814, 816, 818, 820, 822, 824, 826, 828, 830, 832, 834, 836, 838, 840, 842, 844, 846, 848, 850, 852, 854, 856, 858, 860, 862, 864, 866, 868, 870, 872, 874, 876, 878, 880, 882, 884, 886, 888, 890, 892, 894, 896, 898, 900, 902, 904, 906, 908, 910, 912, 914, 916, 918, 920, 922, 924, 926, 928, 930, 932, 934, 936, 938, 940, 942, 944, 946, 948, 950, 952, 954, 956, 958, 960, 962, 964, 966, 968, 970, 972, 974, 976, 978, 980, 982, 984, 986, 988, 990, 992, 994, 996, 998, 999, 1000, 1001, 1002, 1003, 1004, 1005, 1006, 1007, 1008, 1009, 1010, 1011, 1012, 1013, 1014, 1015, 1016, 1017, 1018, 1019, 1020, 1021, 1022, 1023, 1024, 1025, 1026, 1027, 1028, 1029, 1030, 1031, 1032, 1033, 1034, 1035, 1036, 1037, 1038, 1039, 1040, 1041, 1042, 1043, 1044, 1045, 1046, 1047, 1048, 1049, 1050, 1051, 1052, 1053, 1054, 1055, 1056, 1057, 1058, 1059, 1060, 1061, 1062, 1063, 1064, 1065, 1066, 1067, 1068, 1069, 1070, 1071, 1072, 1073, 1074, 1075, 1076, 1077, 1078, 1079, 1080, 1081, 1082, 1083, 1084, 1085, 1086, 1087, 1088, 1089, 1090, 1091, 1092, 1093, 1094, 1095, 1096, 1097, 1098, 1099, 1100, 1101, 1102, 1103, 1104, 1105, 1106, 1107, 1108, 1109, 1110, 1111, 1112, 1113, 1114, 1115, 1116, 1117, 1118, 1119, 1120, 1121, 1122, 1123, 1124, 1125, 1126, 1127, 1128, 1129, 1130, 1131, 1132, 1133, 1134, 1135, 1136, 1137, 1138, 1139, 1140, 1141, 1142, 1143, 1144, 1145, 1146, 1147, 1148, 1149, 1150, 1151, 1152, 1153, 1154, 1155, 1156, 1157, 1158, 1159, 1160, 1161, 1162, 1163, 1164, 1165, 1166, 1167, 1168, 1169, 1170, 1171, 1172, 1173, 1174, 1175, 1176, 1177, 1178, 1179, 1180, 1181, 1182, 1183, 1184, 1185, 1186, 1187, 1188, 1189, 1190, 1191, 1192, 1193, 1194, 1195, 1196, 1197, 1198, 1199, 1200, 1201, 1202, 1203, 1204, 1205, 1206, 1207, 1208, 1209, 1210, 1211, 1212, 1213, 1214, 1215, 1216, 1217, 1218, 1219, 1220, 1221, 1222, 1223, 1224, 1225, 1226, 1227, 1228, 1229, 1230, 1231, 1232, 1233, 1234, 1235, 1236, 1237, 1238, 1239, 12310, 12311, 12312, 12313, 12314, 12315, 12316, 12317, 12318, 12319, 12320, 12321, 12322, 12323, 12324, 12325, 12326, 12327, 12328, 12329, 12330, 12331, 12332, 12333, 12334, 12335, 12336, 12337, 12338, 12339, 12340, 12341, 12342, 12343, 12344, 12345, 12346, 12347, 12348, 12349, 12350, 12351, 12352, 12353, 12354, 12355, 12356, 12357, 12358, 12359, 12360, 12361, 12362, 12363, 12364, 12365, 12366, 12367, 12368, 12369, 12370, 12371, 12372, 12373, 12374, 12375, 12376, 12377, 12378, 12379, 12380, 12381, 12382, 12383, 12384, 12385, 12386, 12387, 12388, 12389, 12390, 12391, 12392, 12393, 12394, 12395, 12396, 12397, 12398, 12399, 123100, 123101, 123102, 123103, 123104, 123105, 123106, 123107, 123108, 123109, 123110, 123111, 123112, 123113, 123114, 123115, 123116, 123117, 123118, 123119, 123120, 123121, 123122, 123123, 123124, 123125, 123126, 123127, 123128, 123129, 123130, 123131, 123132, 123133, 123134, 123135, 123136, 123137, 123138, 123139, 123140, 123141, 123142, 123143, 123144, 123145, 123146, 123147, 123148, 123149, 123150, 123151, 123152, 123153, 123154, 123155, 123156, 123157, 123158, 123159, 123160, 123161, 123162, 123163, 123164, 123165, 123166, 123167, 123168, 123169, 123170, 123171, 123172, 123173, 123174, 123175, 123176, 123177, 123178, 123179, 123180, 123181, 123182, 123183, 123184, 123185, 123186, 123187, 123188, 123189, 123190, 123191, 123192, 123193, 123194, 123195, 123196, 123197, 123198, 123199, 123200, 123201, 123202, 123203, 123204, 123205, 123206, 123207, 123208, 123209, 123210, 123211, 123212, 123213, 123214, 123215, 123216, 123217, 123218, 123219, 123220, 123221, 123222, 123223, 123224, 123225, 123226, 123227, 123228, 123229, 123230, 123231, 123232, 123233, 123234, 123235, 123236, 123237, 123238, 123239, 123240, 123241, 123242, 123243, 123244, 123245, 123246, 123247, 123248, 123249, 123250, 123251, 123252, 123253, 123254, 123255, 123256, 123257, 123258, 123259, 123260, 123261, 123262, 123263, 123264, 123265, 123266, 123267, 123268, 123269, 123270, 123271, 123272, 123273, 123274, 123275, 123276, 123277, 123278, 123279, 123280, 123281, 123282, 123283, 123284, 123285, 123286, 123287, 123288, 123289, 123290, 123291, 123292, 123293, 123294, 123295, 123296, 123297, 123298, 123299, 123300, 123301, 123302, 123303, 123304, 123305, 123306, 123307, 123308, 123309, 123310, 123311, 123312, 123313, 123314, 123315, 123316, 123317, 123318, 123319, 123320, 123321, 123322, 123323, 123324, 123325, 123326, 123327, 123328, 123329, 123330, 123331, 123332, 123333, 123334, 123335, 123336, 123337, 123338, 123339, 123340, 123341, 123342, 123343, 123344, 123345, 123346, 123347, 123348, 123349, 123350, 123351, 123352, 123353, 123354, 123355, 123356, 123357, 123358, 123359, 123360, 123361, 123362, 123363, 123364, 123365, 123366, 123367, 123368, 123369, 123370, 123371, 123372, 123373, 123374, 123375, 123376, 123377, 123378, 123379, 123380, 123381, 123382, 123383, 123384, 123385, 123386, 123387, 123388, 123389, 123390, 123391, 123392, 123393, 123394, 123395, 123396, 123397, 123398, 123399, 123400, 123401, 123402, 123403, 123404, 123405, 123406, 123407, 123408, 123409, 123410, 123411, 123412, 123413, 123414, 123415, 123416, 123417, 123418, 123419, 123420, 123421, 123422, 123423, 123424, 123425, 123426, 123427, 123428, 123429, 123430, 123431, 123432, 123433, 123434, 123435, 123436, 123437, 123438, 123439, 123440, 123441, 123442, 123443, 123444, 123445, 123446, 123447, 123448, 123449, 123450, 123451, 123452, 123453, 123454, 123455, 123456, 123457, 123458, 123459, 123460, 123461, 123462, 123463, 123464, 123465, 123466, 123467, 123468, 123469, 123470, 123471, 123472, 123473, 123474, 123475, 123476, 123477, 123478, 123479, 123480, 123481, 123482, 123483, 123484, 123485, 123486, 123487, 123488, 123489, 123490, 123491, 123492, 123493, 123494, 123495, 123496, 123497, 123498, 123499, 123500, 123501, 123502, 123503, 123504, 123505, 123506, 123507, 123508, 123509, 123510, 123511, 123512, 123513, 123514, 123515, 123516, 123517, 123518, 123519, 123520, 123521, 123522, 123523, 123524, 123525, 123526, 123527, 123528, 123529, 123530, 123531, 123532, 123533, 123534, 123535, 123536, 123537, 123538, 123539, 123540, 123541, 123542, 123543, 123544, 123545, 123546, 123547, 123548, 123549, 123550, 123551, 123552, 123553, 123554, 123555, 123556, 123557, 123558, 123559, 123560, 123561, 123562, 123563, 123564, 123565, 123566, 123567, 123568, 123569, 123570, 123571, 123572, 123573, 123574, 123575, 123576, 123577, 123578, 123579, 123580, 123581, 123582, 123583, 123584, 123585, 123586, 123587, 123588, 123589, 123590, 123591, 123592, 123593, 123594, 123595, 123596, 123597, 123598, 123599, 123600, 123601, 123602, 123603, 123604, 123605, 123606, 123607, 123608, 123609, 123610, 123611, 123612, 123613, 123614, 123615, 123616, 123617, 123618, 123619, 123620, 123621, 123622, 123623, 123624, 123625, 123626, 123627, 123628, 123629, 123630, 123631, 123632, 123633, 123634, 123635, 123636, 123637, 123638, 123639, 123640, 123641, 123642, 123643, 123644, 123645, 123646, 123647, 123648, 123649, 123650, 123651, 123652, 123653, 123654, 123655, 123656, 123657, 123658, 123659, 123660, 123661, 123662, 123663, 123664, 123665, 123666, 123667, 123668, 123669, 123670, 123671, 123672, 123673, 123674, 123675, 123676, 123677, 123678, 123679, 123680, 123681, 123682, 123683, 123684, 123685, 123686, 123687, 123688, 123689, 123690, 123691, 123692, 123693, 123694, 123695, 123696, 123697, 123698, 123699, 123700, 123701, 123702, 123703, 123704, 123705, 123706, 123707, 123708, 123709, 123710, 123711, 123712, 123713, 123714, 123715, 123716, 123717, 123718, 123719, 123720, 123721, 123722, 123723, 123724, 123725, 123726, 123727, 123728, 123729, 123730, 123731, 123732, 123733, 123734, 123735, 123736, 123737, 123738, 123739, 123740, 123741, 123742, 123743, 123744, 123745, 123746, 123747, 123748, 123749, 123750, 123751, 123752, 123753, 123754, 123755, 123756, 123757, 123758, 123759, 123760, 123761, 123762, 123763, 123764, 123765, 123766, 123767, 123768, 123769, 123770, 123771, 123772, 123773, 123774, 123775, 123776, 123777, 123778, 123779, 123780, 123781, 123782, 123783, 123784, 123785, 123786, 123787, 123788, 123789, 123790, 123791, 123792, 123793, 123794, 123795, 123796, 123797, 123798, 123799, 123800, 123801, 123802, 123803, 123804, 123805, 123806, 123807, 123808, 123809, 123810, 123811, 123812, 123813, 123814, 123815, 123816, 123817, 123818, 123819, 123820, 123821, 123822, 123823, 123824, 123825, 123826, 123827, 123828, 123829, 123830, 123831, 123832, 123833, 123834, 123835, 123836, 123837, 123838, 123839, 123840, 123841, 123842, 123843, 123844, 123845, 123846, 123847, 123848, 123849, 123850, 123851, 123852, 123853, 123854, 123855, 123856, 123857, 123858, 123859, 123860, 123861, 123862, 123863, 123864, 123865, 123866, 123867, 123868, 123869, 123870, 123871, 123872, 123873, 123874, 123875, 123876, 123877, 123878, 123879, 123880, 123881, 123882, 123			

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AT	Austria	GB	United Kingdom	MR	Mauritania
AU	Australia	GE	Georgia	MW	Malawi
BB	Barbados	GN	Guinea	NE	Niger
BE	Belgium	GR	Greece	NL	Netherlands
BF	Burkina Faso	HU	Hungary	NO	Norway
BG	Bulgaria	IE	Ireland	NZ	New Zealand
BJ	Benin	IT	Italy	PL	Poland
BR	Brazil	JP	Japan	PT	Portugal
BY	Belarus	KE	Kenya	RO	Romania
CA	Canada	KG	Kyrgyzstan	RU	Russian Federation
CF	Central African Republic	KP	Democratic People's Republic of Korea	SD	Sudan
CG	Congo	KR	Republic of Korea	SE	Sweden
CH	Switzerland	KZ	Kazakhstan	SI	Slovenia
CI	Côte d'Ivoire	LI	Liechtenstein	SK	Slovakia
CM	Cameroon	LK	Sri Lanka	SN	Senegal
CN	China	LU	Luxembourg	TD	Chad
CS	Czechoslovakia	LV	Larvia	TG	Togo
CZ	Czech Republic	MC	Monaco	TJ	Tajikistan
DE	Germany	MD	Republic of Moldova	TT	Trinidad and Tobago
DK	Denmark	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	US	United States of America
FI	Finland	MN	Mongolia	UZ	Uzbekistan
FR	France			VN	Viet Nam
GA	Gabon				

APPARATUS AND METHOD FOR DETERMINING SUBSTANCES  
CONTAINED IN A BODY FLUID

BACKGROUND OF THE INVENTION

The present invention relates generally to an assay system for biological and nonbiological fluids. More particularly, the present invention relates to an apparatus for separating serum or plasma from particulate matter and then optoelectronically evaluating the serum or plasma in order to measure analytes within the serum.

It has long been desirable to utilize devices that can be used for on-site testing of blood products. Particularly important is the analysis of body fluids from humans and animals to diagnose disease, monitor the course of therapy, or determine the presence of illicit drugs. Commonly, the analytical methods used to carry out these objects are performed on blood samples.

Clinical chemists have a preference for working with serum over plasma and plasma over whole blood because of the clarity of the sample matrix and the lack of interfering substances from the solid portion of the blood. In order to facilitate analysis, a separation step must be carried out since the presence of red blood cells, either intact or hemolyzed interferes with the signal generated by the chemical reaction performed by the test.

Conventionally, the separation of blood components has been carried out by placing a blood sample in a centrifuge and centrifuging the sample for ten minutes at approximately 3,000 rpm. The serum obtained from this centrifuging step is then used to carry out the test, thus avoiding interferences from blood solids such as red blood cells.

An embodiment for chemical tests called dry reagent strips was developed first for urinalysis. Thereafter, various efforts to combine dry reagent strip technology in blood testing were started in the early 1950's. Notably, U.S. Patent No. 3,092,465

discloses a reagent in a bibulous carrier with a superimposed semipermeable coating to exclude the chemical and nonchemical interferences from red blood cells. The device, while performing analysis on whole blood, still required additional manipulations by the user, in the form of washing of excess blood after a specified time interval. Additionally, U.S. Patent Nos. 3,552,925 and 3,552,928 disclose the use of salts and amino acids to perform in-situ separation. U.S. Patent No. 4,477,575 discloses the use of a glass fiber matrix.

More recently, membranes have been employed in a variety of devices. These include devices disclosed in the following United States and foreign patents and publications: U.S. Patent Nos. 4,774,192 and 5,166,051; European Published Applications EP 0408222 A1, EP 0408223 A1, EP 0407800 A2 and EP 0388782; and PCT Published Applications Nos. WO 93/22453 and WO 90/10869. The use of the various membranes disclosed in the above patent documents operate on size exclusion principles, and several of these are limited by rates of capillary flow and do not completely eliminate interference from intact or hemolyzed red blood cells. Fresh red blood cells are elastic in nature and may pass through pores smaller than their nominal diameter. Hemolysis may occur on contact with some of the architectural or chemical components of the strips. Consequently, errors may be introduced into the measurement system.

U.S. Patent No. 5,104,619 discloses a disposable diagnostic system comprising a test card having a substantially flat body and a generally cylindrical reagent pad pocket formed in a central area of the flat body. A reagent chemistry pad is disposed in the pocket and a snap fit cover is received in the pocket and arranged over the pad to retain the pad in position. The device size and configuration allows for bar code graphics to be printed on the underneath side of the device. The bar code may contain lot specified data about the reagent chemistry,

and is read by the meter during device insertion. This data may further contain critical parameters for the software algorithm within the meter electronics. U.S. Patent No. 5,139,685 also discloses a separation filter assembly having a snap fit lid.

5 In this patent, glass fibers are utilized and maintained in a compressed state under pressure.

Accordingly, a need exists for an integrated system for assaying analytes and whole blood samples which are not affected by the chemical or physical interferences normally caused by red 10 blood cells and other portions of whole blood.

#### SUMMARY OF THE INVENTION

The present invention, in one form thereof, comprises a dry solid phase diagnostic test strip and system for the chemical, enzymatic, and/or immunological analysis of whole blood 15 analytes, comprising a reflectance photometer, a solid support strip, a porous detection zone member, a permeable spreading layer, an overlay sample receiving membrane containing an agent for the exclusion of intact red blood cells and a strip-receiving platform for positioning the strip inside the 20 reflectance photometer. The detection area membrane may contain chemical, enzymatic, and/or immunological reagents that generate specific signals in the presence of a target analyte. The agent, in contact with the overlay membrane, prevents passage and hemolysis of red blood cells while facilitating rapid 25 transport and reaction of the plasma or serum portion of introduced whole blood samples.

In addition, the present invention, in one form thereof, comprises a reflectance photometer which utilizes test strips that are color coded for test differentiation. For example, a 30 blue strip may indicate a glucose test, whereas a red strip may indicate a cholesterol test. These colors are then divided into shades such as 64 shades of blue equal to 64 lot numbers of glucose strips. The photometer includes a separate optical read

head that determines the color and shade of the base of the test strip device as the strip is inserted into the photometric instrument. The shade is converted into a lot number ranging from 1 to 64. The instrument also has a memory module (preferably an electrically erasable programmable read-only memory) that has a corresponding lot number to the shade of the strip to ensure lot number verification. The instrument then compares the inserted memory module programmed lot number to ensure that it is the same lot number as the test strip. If the strip lot number does not match the memory module lot number, the test is not performed, and the user is instructed to insert the correct memory module.

The lot number verification allows for the automated coding of lot numbers so that the user does not need to enter a lot code for each vial of strips. This prevents the running of the incorrect, old, or expired lot number tests in the instrument.

The "plug-in memory" of the module includes the lot number of the test strip, the expiration date, and the performance criteria for the actual strip measurement. The performance criteria include the wavelength, measurement algorithm, and unreacted density qualifications necessary for a valid test result.

The optoelectronic measurements of the chemistry test reaction on and in a surface enhances the dynamic range of the dry phase test. Algorithms that read at different wavelengths at different times in the chemistry reaction can extend the dynamic range of the test system. This is particularly applicable when using multiple chromophores in a single measurement system. The early portion of a chemistry could be read at the peak wavelength of a reaction, while the later portion or darker or more dense portion of color development could be read at a wavelength not near the peak of the color development. In addition, different chromophores may respond in

a more linear manner in different portions of the dynamic range of the chemistry. Manipulation of these two data points can significantly increase the dynamic range (in mg/dl) of a chemistry reaction.

5        The optoelectronic measurement of the chemistry test reaction on and in a surface reduces error due to orientation of the surface to be read to the instrument. Multiple wavelengths and different angles are used to correct possible problems in positioning the strip in the instrument. If the detector is at 10 "0" angle and the emitters of the same or different wavelengths are at different angles (e.g., one at 40° and one at 50°), the tilting of a surface will positively contribute to one reading while it will contribute in a negative manner to the other reading thus it is able to cancel the error presented by the 15 angle presentation of the surface. These same measurement methods can be used to eliminate interferences from substances such as bilirubin and others.

20       The optoelectronic measurements of the chemistry test reaction on and in the surface enhance the stability of timed and untimed dry phase chemistry reaction. Algorithms are used to determine the "end point" of a chemistry. In other words, measurements can be done at similar or dissimilar wavelengths to predict the stable portion or end point of a chemistry. If kinetic measurements are made, the kinetic readings can be 25 subjected to an algorithm to determine that the rate is slow enough to declare the extrapolate chemistry is at an end or completion. When known standards are run and predicted by this pseudo-endpoint, the same measuring criteria can be applied to unknowns to determine the "endpoint" of the test reaction.

30       The use of colored or shaded visual indicators in the instrument enhance the interpretation of test results. A colored bar graph is used to aid the user in knowing when the user test results are in a normal or safe range. Out of range

colors such as orange for caution and red for danger are used when results are outside the green "safe" range. This is particularly useful to new testers who are not familiar with the number scale of the different test results. A voice module can 5 also be used to warn the user of unsafe results or operation of the instrument system to make the system usable by the visually impaired by providing, for example, a sound beep for each unit of glucose during a glucose test.

#### BRIEF DESCRIPTION OF THE DRAWINGS

10 Fig. 1 is a perspective view of the reflectance photometer in accordance with an embodiment of the present invention;

Fig. 2 is an exploded perspective view of the plastic test strip of present invention in its unlocked position;

15 Fig. 3 is a perspective view of the plastic strip of Fig. 2 in its locked position; and

Fig. 4 is a sectional view of the plastic strip;

Fig. 5 is a block diagram schematic of one embodiment of the reflectance photometer of the present invention;

20 Fig. 6 is a graph plotting sample size, elapsed test time and percentage of reflectance illustrating how endpoint determinations may be utilized to speed chemistry measurement.

#### DESCRIPTION OF THE PREFERRED EMBODIMENT

In accordance with the embodiment of the present invention, the diagnostic chemistry measurement device 10 for dry solid 25 phase chemical, enzymatic, immunological assay of whole blood or sera analytes is made up of an injection molded carrier test strip 20 in which several porous and nonporous materials containing chemicals and reactants are contained for the purpose of generating a detectable signal in the presence of certain 30 analytes. The test strip 12 is inserted into a reflectance photometer. The reaction material layer on the test strip 12 is held in intimate noncompressed contact with a whole blood separation layer in the absence of adhesives for the purpose of

providing a liquid sample free of red blood cells to the reaction layer or layers.

Holder

The holder test strip 12 of this invention acts as holder for the different layers of the test reaction system. It provides a convenient handle as well as a mechanism for placing test strip 12 into an instrument 10 for the reading of the density changes of the reaction layers. As shown in Fig. 2 test strip 12 includes an elongate body 16 preferably formed by injection molding. Elongated body 16 includes a first end portion 18 and a second end portion 20. A hinged portion 22 is located between first and second end portions 18 and 20 so that first end 18 is foldable over elongated body 16 into contact with second end 20.

As shown in Fig. 2 first end portion 18 includes a opening 24 while second end portion 20 includes a complementary spaced opening 26. When first end portion 18 is folded over body 16, each opening 24 and 26 are aligned. In its folded position as shown in Fig. 3 opening 24 in test strip 12 defines an area for depositing a body fluid sample while opening 26 defines an area in which optoelectronic measurements of chemistry test reactions are conducted.

Test strip 12 further includes an adhesiveless carrier layer 14 formed from, for example, three particular layers. In a standard diagnostic test strip, carrier layer 14 may include a disbursement layer 28, formed of for example woven materials such as polyester or cotton, for rapid and even disbursement of body fluid along carrier layer 14. Beneath that may be included a separating layer 30 constructed of known materials such as shown in Table IX infra, that, when exposed to a sample liquid, may separate analyte and analyte disrupting elements such as red blood cells from whole blood. This action would permit the serum analytes to pass through separating layer 30 while

5 preventing red blood cells or other analyte disrupting elements from passing through. The last layer shown in Fig. 2 is that of the test reaction membrane 32 on which the dry chemicals and reactants are contained for generating a visible signal in the presence of serum analytes. Molded carrier body 16 serves as a support for the reacting and nonreacting layers 28, 30 and 32 which may be formed from papers, membranes and deles materials.

10 The test strip holder 12 positions the different layer materials 28, 31, 32 within the holder the correct x, y, and z axis positions. Carrier layer 14 made up, for example, the disbursement separating and test reaction layers 28, 30 and 32 are held in noncompressed adhesiveless locations by first end portion 18 folding over to second end portion 20. This may be accomplished in a number of different ways. The preferred way 15 of noncompressingly holding carrier layer is of an upstanding annular rim 34 may help locate the carrier layer 14 within test strip 12. Additionally, small upstanding protuberances 36 along second end portion 20, radially located away from opening 26 prevent movement of carrier layer 14. The purpose of both 20 annular rim 34 and small upstanding protuberances 36 is to hold the layers of carrier layer 14 without compression between opening 24 and opening 26, thereby preventing pooling of any sample within carrier layer 14. This consideration of noncompression of the carrier layer 14 is of greater importance 25 when larger numbers of layers are utilized. The positioning of a carrier layer 14 without adhesives or compression allows for efficient transport of sample and reactants contained in the system and test strip 12. Annular rim 34 or alternatively other areas of test strip 12 may include sawtooth protrusions to 30 increase flow rate through carrier layer 14.

Test strip 12 includes a locking mechanism to prevent any unlocking of front end portion 18 from its folded position over elongated body 16. As shown in Fig. 2, one type of locking

mechanism may include a plurality of upwardly extending tabs or projections 38 that interfit or lock into corresponding openings 40 in first end portion 18. When first end portion 18 is folded to second end portion 20, lock projections 38 will interfit and snap lock within openings 40. Other types of one way locking mechanisms may also be used, such as snap rivets.

More than one test reaction system can be housed in a test strip 12. A second set of holes 24, and 26 may be included in test strip 12 so that two tests may be run at once.

10 The described holding mechanism allows for the rapid separation of whole blood into its liquid and solid components. It also allows sample volumes as low 2.0 microliters to be used in dry phase chemistry reactions. Test strip 12 allows the use of several reaction and non-reaction layers. A typical holder 15 could contain from 1 to 8 layers of material with thicknesses from approximately 0.002 inches to 0.007 inches, for example.

Chemicals and materials are employed to allow for the treatment of samples such as whole blood, which will allow the whole blood sample to be separated without disrupting the red 20 blood cells while rapidly moving the liquid portion of the whole blood sample to one or more reaction sites in the holder, normally on a test reaction membrane 32. These chemicals can be composed of polymeric and nonpolymeric substances that are dried onto one or more surfaces of the materials contained in the 25 device holder. Additionally, light metal salts of elements such as Potassium, Lithium, Sodium, and Calcium may be utilized to treat red blood cells before and during the separation process. The materials which may be used in the holder for treatment by or containment of these chemicals can be composed of woven, 30 nonwoven, napped, or flocked materials.

Analytes

A wide variety of analytes can be determined by using the disclosed apparatus. Examples are given in tables I and II, *infra*.

- 5        Further, given the small size and robust nature of the reagent strips and instrumentation, analyses need not be limited to traditional clinical laboratory settings. The device of the present invention is also simple enough to be used by people with minimal or no chemical or medical technology training.
- 10      This advantage allows use at home, or by mobile health care delivery people. Examples of this are diabetics that must monitor themselves for glucose and ketone bodies, patients on home dialysis who would benefit by monitoring of urea nitrogen and people endeavoring to lower their cholesterol levels.
- 15      Further, by combining several different reagents on a single support, a panel of tests may be done. Examples of this would be a liver panel consisting of ALT, AST, Alkaline Phosphates. A diabetic panel might consist of glucose, beta hydroxybutyrate and glycated hemoglobin. A coagulation panel
- 20      might consist of Prothrombin time, ACTT, and ACT.

## FAMILIES OF ANALYTES BY STRUCTURE

Table I

Family	Examples
Carbohydrate	glucose, lactose, galactose
Nitrogen Molety	urea nitrogen, creatinine, uric acid
Lipid	cholesterol, triglycerides, LDL, HDL
Enzyme	ALT, AST, Alkaline Phosphatase, CPK, CK-MB
Hormone	HCG, LH
Therapeutic Drugs	theophylline
Drugs of abuse	cocaine, marijuana, barbiturates, salicylates
Electrolyte	Na <sup>+</sup> , K <sup>+</sup> , Cl <sup>-</sup> , Li <sup>+</sup> , CO <sup>2</sup>
Nucleic Acids	infectious disease, forensic applications, genetic disorders

## FAMILIES OF ANALYTES BY DISEASE

Table II

Disease	Examples
Diabetes	glucose, beta hydroxybutyrate, hemoglobin A <sub>1c</sub>
Liver problems	ALT, AST, bilirubin
Acidosis/Alkalosis	PO <sub>2</sub> , PCO <sub>2</sub> , pH
Hypertension	Na <sup>+</sup> , K <sup>+</sup>
Nutritional status	Ca <sup>++</sup> , Mg <sup>++</sup> , Zn <sup>++</sup> , trace minerals

Examples

The following illustrative examples teach various combinations of buffers, dyes, stabilizers and other reactive and functional components which may be combined by a person having ordinary skill in the art into the system test reaction areas.

Table IX gives various types of dyes and indicators used in diagnostic reagents.

## Example #1

## 10 Glucose measuring system

TABLE III

Ingredient	Function	Amount	Available from
Glucose Oxidase	reactant	25,000	Sigma Chemicals, St. Louis, MO
Peroxidase	reactant	75,000	Sigma Chemicals, St. Louis, MO
Silwet 7500	surfactant	0.10ml	Dow-Corning, Midland, MI
PVP K 30	enzyme stabilizer	0.50gms	ISP, Linden, NJ
Citric Acid	Buffer system	1.25gms	Aldrich Chemical, Milwaukee, WI
Sodium citrate	Buffer system	0.10ml	Dow-Corning, Midland, MI
DOW 1520	antifoam	1.00gms	Aldrich Chemical, Milwaukee, WI

4 AAP	chromophore	0.25gms	Aldrich Chemical, Milwaukee, WI
3,5 DCHBS	chromophore	0.25mgs	Boehringer Mannheim
Distilled H <sub>2</sub> O	solvent	QS to 100ml	

Preparation: Approximately 50ml of distilled H<sub>2</sub>O was placed in a beaker on a stirring plate. A magnetic bar was added and the ingredients added sequentially after the previous gradient was dissolved and dispersed. After all ingredients were added the volume was adjusted to 100ml of distilled H<sub>2</sub>O.

5

#### Example 2

##### Triglycerides measuring system



**TABLE IV**

Ingredient	Function	Amount	Available from
cholesterol esterase	reactant	15,000 units	Shinko-American, N.Y., N.Y.
glycerol kinase	reactant	5,000 units	Shinko-American, N.Y., N.Y.
glycerophosphate oxidase	reactant	5,000 units	Shinko-American, N.Y., N.Y..
peroxidase	reactant	5,000 units	Shinko-American, N.Y., N.Y.
4 AAP	chromogen	1.00gm	Aldrich
3, 5 DCHBS	chromogen	0.25gm	Boehringer Mannheim
MES	buffer	2.50gm	Research Organics
PVP K30	stabilizer	0.50gm	ISP
glucose	filler	2.50gm	Sigma
triton X-100	surfactant	0.10gm	Boehringer Mannheim
Distilled H <sub>2</sub> O	solvent	QS to 100ml	

Preparation: Same as example #1

Example 3

Cholesterol measuring system (all amounts approximate)

TABLE V

Ingredient	Function	Amount	Available from
cholesterol Oxidase	reactant	10,000	Shinko-American, N.Y., N.Y.
cholesterol esterase	reactant	7,000	Shinko-American, N.Y., N.Y.
sodium phosphate 0.5 M pH 7.0	buffer	750 ml	Dow-Corning,
B.S.A.	surfactant	15 gm	Aldrich Chemical,
peroxidase	reactant	170,000	Shinko-American,
DOSS	surfactant	7.0 gms	Boehringer Mannheim
sucrose	stabilizer	1.0 gms	Sigma Chemicals,
TMB	chromogen	10.0 gms	Aldrich Chemical,
Distilled H <sub>2</sub> O	solvent	QS to 100 ml	

Preparation: same as example #1

Alternatively, the chromogen may be prepared in an organic solvent matrix and treated as a first or 2nd application to the membrane or paper.

5

TABLE VI

Ingredient	Function	Amount	Available from
Acetone/methanol 1:1	solvent	100ml	Aldrich
Tetramethyl benzidine	solvent chromogen	1.00gm	Biosynth Inc., Chicago, IL

Example 4

Blood Urea Nitrogen Measuring System

**TABLE VII**

Ingredient	Function
Urease	reactant
H <sub>2</sub> O	solvent
Bromothymol blue	chromogen
PVP K90	film former
Fructose	filler

Preparation: Same as experiment #1.

**Table VIII****Types of Indicators****Chromogenic substrate**

5                    Redox

Leuco dyes

Oxidative couplers

Benzidene Derivatives

Fluorescent labels

10                Dye releasing system

**Table IX****Separation mechanisms used in dry reagents**

Chemical	Physical	Mechanical
Dextran	hydrophilic polymers	centrifuge
sugars	porous latex films	filters
lectin	polymer & swelling agent	filters & pressure
amino acids	membranes	membranes & differential pressure
PEG/polyacrylate	microfiber cloth	wedge shape
thrombin	napped cloth	
gels	sintered porous matrix	
coagulants	density gradient	
agglutinating agents	glass fibers	

amine polymers	hollow fibers	
trivalent cations	membrane	

Spectrophotometer

The present invention also includes use of a spectrophotometric device 10 for determining the density of the color reaction on and in the membrane surface of the test reaction layer 32 within test strip 12. Photometric device 10 as shown in Fig. 1 includes a hand-held housing 50 for containing electronic control circuitry for operating the aforementioned tests. In the embodiment shown in Fig. 1, a test strip holding region 52 is located above three light detectors or sensors 54 each disposed within a port 56. During test operation, a test strip 12 is inserted into holding region 52 so that test strip openings 26 are located adjacent ports 56. Light sensors may take a reading from light reflected from the exposed test reaction membrane layer 32 or from test strip 12 itself to determine its color.

Housing 50 further includes a specialized display device, such as a liquid crystal display 58. Display 58 is utilized for relating test results and other information to the user. In particular, a color scale 60 is used to facilitate interpretation of test results operating concurrently with digital display segments 62. Additional display segments on display 58 include a test wait indicator segment 64 to inform the user to wait while device 10 is performing the selected tests, and a test name segment 66 which the unit determined from the type of test strip 12 inserted.

Color scale 60 may easily be constructed by a plurality of shaded or colored segments arranged adjacent each other to form a bar graph like indicator. Electrically controllable segments 68 are oriented over the color or shaded segments so that when segments 68 are activated segments 68 become dark, preventing

certain colored or shaded segments 60 from being visualized or viewed. Segments 68 that are not activated permit the underlying colored or shaded segments of color scale 60 to be visualized. In this way it is possible for an electronic control to permit only a single colored or shaded segment to be viewed thereby communicating test results.

5 A possible result range spectrum for color scale indication segments may include particular colors with particular test result meanings such as:

- 10 very high result danger, RED
- high result danger, RED
- high result caution, YELLOW
- high result caution, YELLOW
- high normal result, GREEN
- 15 normal result, GREEN
- normal result, GREEN
- low normal result, GREEN
- low result caution, YELLOW
- low result caution, YELLOW
- 20 very low result danger, RED

Color scale 60 permits an unsophisticated user to instantly visually determine, in one embodiment, if a test result is normal (a green segment visualized), slightly abnormal (a yellow segment visualized) or dangerous high or low result (a red segment visualized). Alternatively, if a color liquid crystal display is utilized, the electronic control for test unit 10 may directly indicate a colored segment, rather than covering all 25 but one colored segment.

30 A suitable instrument, such as a diffuse reflectance spectrophotometer 10 with appropriate software, can be made to automatically read reflectance at certain points in time, calculate the rate of reflectance change, and by using calibration factors and software, output the level of analyte in

the fluid tested. The electronic control mechanism of photometric unit 10 is shown in schematic form in Fig. 5. One or more light sources 70, for example high intensity light emitting diodes (LED) are disposed in housing 50 to illuminate test strip 12 as shown by arrows 72. A light detector or sensor 54, for example a photo transistor, is able to take a reading of light reflected either from the surface of test strip 12 or from its associated test reaction membrane 32. Light source 70 and light sensor 54 can be adapted to generate or respond to particular wavelengths of light.

Sensor 70 transmits a signal to an amplifier 74 as is known in the art. One type of amplifier available for use is, for example, a linear integrated circuit which converts the phototransistor current to a voltage signal.

Appropriate electronic circuitry is utilized to take the output of amplifier 74, normally a sample and hold unit 76, and transfer the signal to an analog-to-digital converter 78. Analog-to-digital converter takes the analog voltage output from the sample and hold unit 76 and converts it to, for example a 16 bit binary digital number upon command of a microprocessor/microcontroller unit 80.

Preferably an electronic microprocessor/microcontroller 80 utilizing digital integrated circuitry is used to time selected tests, read signals, and together with associated programs and data memory 82, calculate and store reflectivity values and calculate analyte levels from the stored data.

Additional information for particular tests may be stored in a removable EEPROM unit 84 operably connected to microprocessor/microcontroller 80. EEPROM unit 84 is an interchangeable plug-in memory module containing measurement parameters, software, calibration data, and reagent recognition data for particular test strips 12. Additionally, EEPROM unit 84 contains the shelf life data and identity verification

information for particular production runs or lots of test strips 12.

Automated lot coding is done by the color coding the plastic material used to make the test strip holder 12. The 5 color used in test strip holder 12 preferably has 16 different densities that can be distinguished by at least one of the wavelengths used in the optical sensor head 54 of instrument 10. For instance the dynamic range of the  $\pm$  reflectances of the strip holder color could be as follows to determine the 10 different shades of color density:

	$\pm$ Reflectance Green LED	$\pm$ Reflectance Red LED	Lot #
15	70		1
	65		2
	55		3
	50		4
	45		5
	40		6
20	35		7
	30		8
		70	9
		65	10
		60	11
		55	12
25		50	13
		45	14
		40	15
		35	16

As the strips 12 are inserted into device 10, the instrument 30 detects a change in the measurement area. This change indicates that a strip 12 has been inserted into the instrument 10. As the instrument detects the insertion of a test strip 12, it reads the densities of at least one of the LED's and calculates the lot number by the above table. Instrument 10 then goes to 35 the EEPROM port connected to microprocessor / microcontroller 80 which has an EEPROM unit 84 inserted. Instrument 10 checks to see that the EEPROM preselected lot number is the same as lot number of test strip 12 that had been inserted into the instrument. If the lot numbers are the same for test strip 12 40 and EEPROM 84, the instrument downloads the information contained in the EEPROM and proceeds with the test analysis.

The instrument 10 reads the density of the unreacted strip to assure quality of the strip before the test is initiated, if quality is passed then the instrument instructs the user to apply a sample.

5 A sample is then applied and instrument 10 begins a measurement cycle to ensure that the proper amount of sample was applied to the test strip. When the instrument has determined that enough sample has been applied, it then goes into another cycle to measure the end of the chemistry reaction. When the 10 end of the chemistry reaction has occurred, then the instrument measures the final density and compares it to a measurement algorithm stored in EEPROM unit 84. This measurement algorithm then determines the concentration of the test to be measured by comparing the measured density (darkness) of the color formed 15 and comparing this density number to a table of values through the use of an algorithm stored in the EEPROM unit 84.

After a particular test strip is selected and placed in the unit, a sample, normally a whole blood sample from a fingertip or from a pipiter tip (which could have gotten its sample from a 20 tube of blood as in a laboratory type situation) is applied to the sample application spot, opening 24, on test strip 12. A dispersemant layer 28 causes the sample to quickly spread over the entire area of carrier layer 14. The separation layer 30 of the test strip spot is allowed to separate out the solids (red 25 blood cells and other analyte disrupting elements) from the liquid (plasma or sera or other analyte containing portion). The separated fluid, i.e the plasma, sera, or other analyte containing portions, moves to the test reaction membrane layer 32 below the separation membrane 30. The above fluid migration 30 causes the reactants (analytes such as glucose) in the sample to come into contact with the reactants in test reaction membrane layer 32.

Analytes/fluid contacts reagent layer reaction 32 and initiates an appearance or disappearance of color, depending on its particular reaction. The above presentation of analyte to the reaction layer 32 causes the desired reaction to occur.

5 This reaction causes a color change that can be detected both visually and by the instrument. The color change is then converted into a digital result on the instrument LCD as described above. A comparison color chart can be used to visually determine a reaction quantity scale as in litmus paper.

10 Instrument 10 can use different wavelengths at different density portions of the reaction to maximize the dynamic range of the chemistry and the limits of the instrument at a particular wavelength.

15 The "end-point" of the reaction is defined as a point where there appears to be no change or a very small change in density. That is, the chemistry changes color proportional to the concentration of the reactants that has come into contact with the reactants materials in the test pad (membrane). This small amount of change can be a change per time period. An example  
20 would be as per the graph in Fig. 6. Detailed information used to generate this graph is that the changes per 5 second time period during the beginning of the test reaction would be greater than 5% reflectance per 5 second time period. When this change is less than 1% reflectance per time period it can be  
25 said that the reaction is complete or at an endpoint. The instrument stores this percentage reflectance at this time and uses as above to determine the concentration of the analyte tested for in the test strip.

30 The Kubelka-Monk equation of  $K/S = (1 - \text{reflectance})^2$  divided by  $(2 \times \text{reflectance})$  can be used to linearize the percentage reflectance values. This linearization simplifies the algorithm necessary to calculate results. This pseudo endpoint chemistry allows a more stable read time, which in turn allows for a more

reproducible answer. Pseudo endpoints also permit a more rapid assay to be performed. Certain other glucose monitoring systems incorporate pre-determined timing circuit. This pseudo endpoint allows for a different method to be used in measuring chemistry reactions, provided one can determine the endpoint of the chemistry by a method other than timing.

Multiple wavelengths are used to enhance the dynamic range of a chemistry. This is particularly useful when one uses a multiple chromophore indicator system as do some of the above mentioned chemistries. Early portions or low concentrations of a test such as glucose can use a broad range indicator such as TMB to increase sensitivity in the low to mid range of the chemistry. When the test concentration is higher or the reaction faster, a different chromophore is focused upon to determine more dynamic range than the previous chromophore. This allows one to expand the dynamic range by two different methods.

One can also use wavelengths on the peak for more dynamic range and wavelengths off the "peak" absorbance of the test system to enhance or reduce dynamic range and also to enhance or reduce the "pseudo endpoint" algorithms. Manipulation of these four factors, chromophore A, chromophore B, wavelength 1 and wavelength 2 can allow one to better define the "pseudo endpoint" algorithm and also allow one to optimize the dynamic range of the chemistry which in turn allows for increased sensitivity throughout the chemistry reaction range with greater precision.

Multiple wavelengths can also be used with different angles of emission to correct possible problems in positioning the strip in the instrument. If the detector is at "0" angle and the emitters of the same or different wavelengths are at different angles (one at 40° and one at 50°) the tilting of a surface will positively contribute to one reading while the

other contributes in a negative manner thus cancelling the error presented by the angle presentation of the surface. These same measurements methods can be used to eliminate interferences from substances such as bilirubin and others. When the angle of 5 light incidence is increased from improper positioning of a chemistry read surface to the instrument optics, errors of both gloss and angularity are introduced into the measuring system and can give false low readings.

Examples

- 10 1. Indicators and chromogens advantageously used in combination
  - a. wide range pH test.  
Bromothymol blue and methyl red covers pH range of 5 through 9
  - b. 4 amino antipyrine + 3,5 dichlorohydroxybenzene sulfonate (4AAP+3,5 DCHBS)
  - c. TMB+Chromotropic acid
  - d. Syringaldazine + Vanillin Azine
- 20 2. Color coding for test and lot identification
  - a. blues, 16 different shades (density)
  - b. reds, 16 different shades (density)
  - c. greens, 16 different shades (density)
  - d. yellows, 16 different shades (density)
  - e. oranges, 16 different shades (density)
  - f. browns, 16 different shades (density)
  - 25 g. magentas, 16 different shades (density)
  - h. light blues, 16 different shades (density)
  - i. light reds, 16 different shades (density)
  - j. light greens, 16 different shades (density)
  - k. light browns, 16 different shades (density)
  - 30 l. light magentas, 16 different shades (density)
  - m. cyan, 16 different shades (density)

n. light cyan, 16 different shades (density)

It will be appreciated that the foregoing is presented by way of illustration only, and not by way of any limitation, and that various alternatives and modifications may be made to the 5 illustrated embodiment without departing from the spirit and scope of the invention.

## WHAT IS CLAIMED IS:

1. A diagnostic test strip for use in an analyzer for measuring analyte in a sample, said test strip comprising:
  - 5 an elongate body including a first end, a second end, and a hinged portion between said first and second ends, said first end being foldable over said body, said first end and said second end each having an opening that are aligned when said first end is folded; and
  - 10 an adhesiveless carrier layer disposed without compression between said first end and said body whereby sample communicated to said carrier layer is prevented from pooling within said carrier layer.
2. The test strip of Claim 1 in which said carrier layer includes a separating layer that when exposed to a whole blood sample excludes red blood cells from passing therethrough while allowing the liquid portion of the whole blood sample to pass therethrough.
3. The test strip of Claim 1 having a carrier layer utilizing samples in the range of 2.0ul to 10.0ul to generate a reaction to accurately test a selected analyte.
4. The test strip of Claim 1 in which said carrier layer includes:
  - 5 a separating layer that when exposed to a sample liquid having both analytes and analyte disrupting elements said separating layer excludes said analyte disrupting elements from passing therethrough while allowing the analyte portion of the sample liquid to pass therethrough; and
  - 10 a test reaction membrane adjacent said separating layer that creates a gradient color dependant on the concentration of selected analytes in the analyte portion that had passed through said separating layer.
5. The test strip of Claim 1 in which said carrier layer

includes a spreading layer over said separating layer to cause sample to substantially evenly enter said separating layer.

6. The test strip of Claim 1 in which said carrier layer may test more than one analyte at one time.

7. The test strip of Claim 1 in which said body includes a locking means lock together said first end and said body.

8. A diagnostic test strip for use in an analyzer for measuring analyte in a sample, said test strip comprising:

an elongate body including a first end, a second end, and a hinged portion between said first and second ends, said first end being foldable over said body, said first end and said second end each having an opening that are aligned when said first end is folded;

an adhesiveless carrier layer disposed between said first end and said body of said carrier layer whereby the accuracy of measured analytes is increased; and

non-compressive means holding said carrier layer in place between said first end and said body whereby pooling of sample within said carrier is prevented.

9. The test strip of Claim 8 in which said non-compressive means comprise protrusions to locate said carrier layer in place whereby said carrier layer is maintained in known locations along the X, Y and Z axis.

10. The test strip of Claim 8 in which said non-compressive means comprise sawtooth protrusions to locate said carrier layer in place whereby said carrier layer is maintained in known locations along the X, Y and Z axis.

11. The test strip of Claim 8 in which said separation layer is treated with light metal salts to reduce red blood cells in the sample.

12. The test strip of Claim 8 in which said first end folds over said carrier layer and locks to said body causing the layers of said carrier layer to be in adjacent contact without

5        adhesives or compression whereby efficient separation of red blood cells from plasma in whole blood samples.

13. The test strip of Claim 1 further comprising a tab and an opening one of which on said first end the other on said second end so that when said first end folds into contact with said body, said tab interfits with said opening to lock said 5 first end with said body.

14. A chemistry measurement system comprising:  
a test instrument with a light source and light sensor;

5        diagnostic test strip for use in analyzing a sample, an elongate body including a first end, a second end, and a hinged portion between said first and second ends, said first end being foldable over said body, said first end and said second end each having an opening that are aligned when said first end is folded, said test strip having an adhesiveless 10 carrier layer disposed without compression between said first end and said body, said opening adapted to receive said sample;

an electronic control for computing particular test results on light incident on said light sensor that was reflected from said test strip; and

15        display means for displaying said test results.

15. The measurement system of Claim 14 in which a plurality of test strips are utilized for particular chemical tests, said test strips color coded for identification of which said particular chemical test said test strip is operable.

16. The measurement system of Claim 15 in which a plurality of test strips are utilized for particular chemical tests, said test strips color coded for identification of said particular chemical test for which said test strip is operable 5 and lot designator, said light sensor sensing said color of said test strip when said strip is inserted into said test instrument and sending a coded signal to said electronic control, said

- 10        electronic control determining the type of test and lot designator of inserted test strip and denying test operation when said lot designator is not within predefined limits.
17. A chemistry measurement system of Claim 16 in which said system tests for more than one analyte with a single test strip.
18. A chemistry measurement system comprising:  
a test instrument with a light source and light sensor;  
a color coded diagnostic test strip for use in  
5        analyzing a sample, said strip having an elongate body including a first end, a second end, and a hinged portion between said first and second ends, said first end being foldable over said body, said first end and said second end each having an opening that are aligned when said first end is folded, said test strip  
10        having an carrier layer disposed without compression between said first end and said body, said first end opening adapted to receive said sample;
- 15        an electronic control for computing particular test results on light incident on said light sensor that was reflected from said test strip, said control determining from said color of said test strip if said test strip is from a particular production lot, said control operating said test operation only if said test strip is from a preselected production lot; and
- 20        display means for displaying said test results if said test is conducted.
19. A chemistry measurement system comprising:  
a test instrument with a light source and light sensor, said light source emitting light at multiple angles and multiple wavelengths;
- 5        a diagnostic test strip for use in analyzing a sample, said test strip comprising an elongate body including a first

end, a second end, and a hinged portion between said first and second ends, said first end being foldable over said body, said test strip having an adhesiveless carrier layer disposed without 10 compression between said first end and said body;

an electronic control for computing particular test results on light incident on said light sensor that was reflected from said test strip; and

display means for displaying said test results.

20. The chemistry measurement system of Claim 19 in which said electronic control includes a removable erasable programmable read only memory unit containing lot number data and expiration data for particularly chemical tests.

21. A liquid crystal display matrix for a hand-held chemistry measuring system, said display matrix comprising:

a display screen having a plurality of shaded segments arranged adjacent each other;

5 a plurality of controllable segments disposed within said screen oriented over said shaded segments, said controllable segments preventing visualization of said shaded segments when activated and permitting visualization of said shaded segments when deactivated whereby results from said 10 chemistry measuring system are communicated by visualizing selected shaded segments.

22. The liquid crystal display of Claim 21 in which said shaded segments are aligned in a line forming a substantial bar graph indicator.

23. The liquid crystal display of Claim 21 in which said shaded segments are colored to indicate selected results from said blood chemistry measuring system.

24. A method of testing analyte in a sample comprising the steps of :

providing a chemistry measurement system having a light sensor, test display, colored diagnostic test strips for

- 5        use in analyzing a sample, each said test strip having an adhesiveless carrier layer attached without compression to said test strip;
- selecting a test strip of a predetermined color for a predetermined test;
- 10        depositing a test sample on the carrier layer of the selected test strip;
- inserting said test strip into said measurement system;
- operating said measurement system to conduct the predetermined test; and
- 15        displaying test results on said test display.

1 / 5

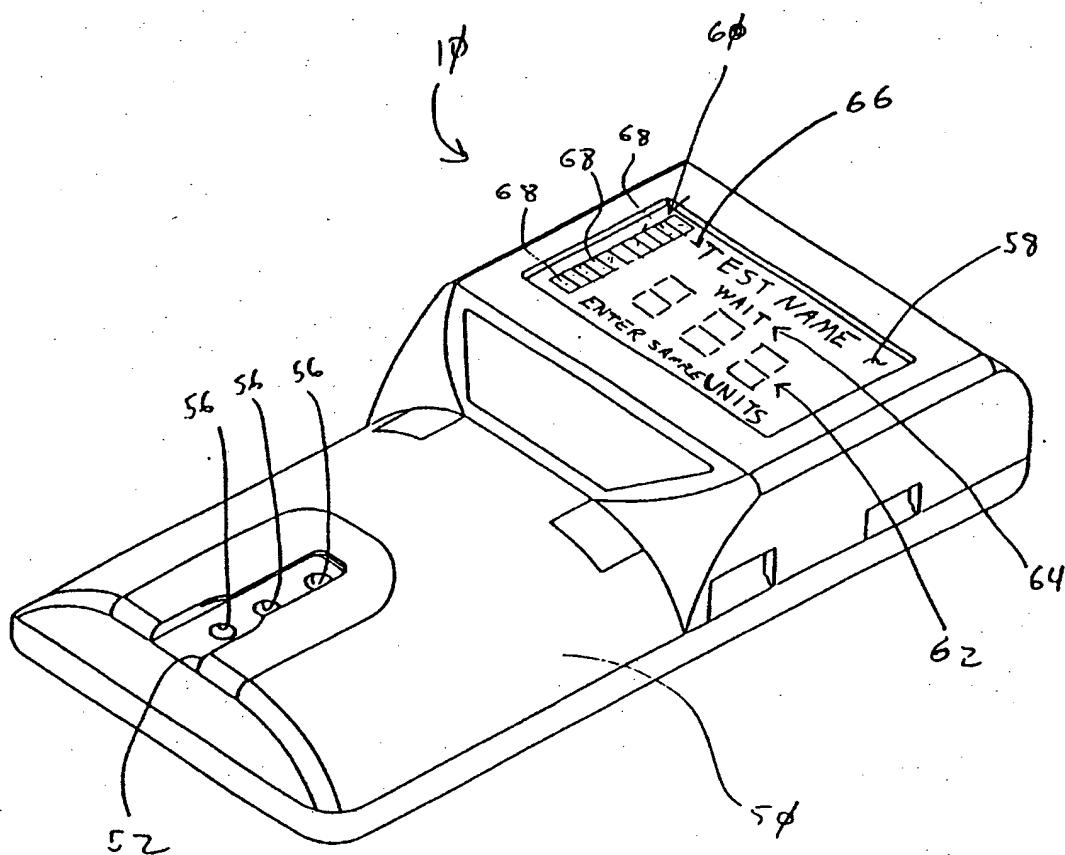


Fig 1

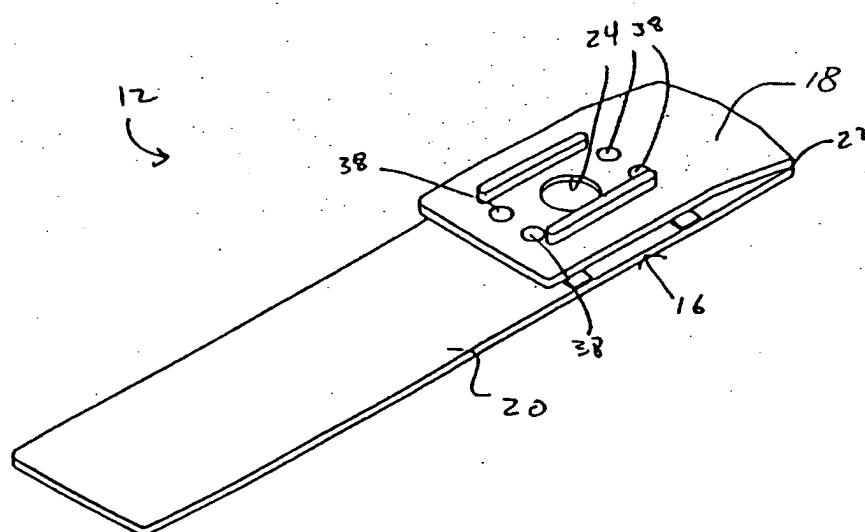
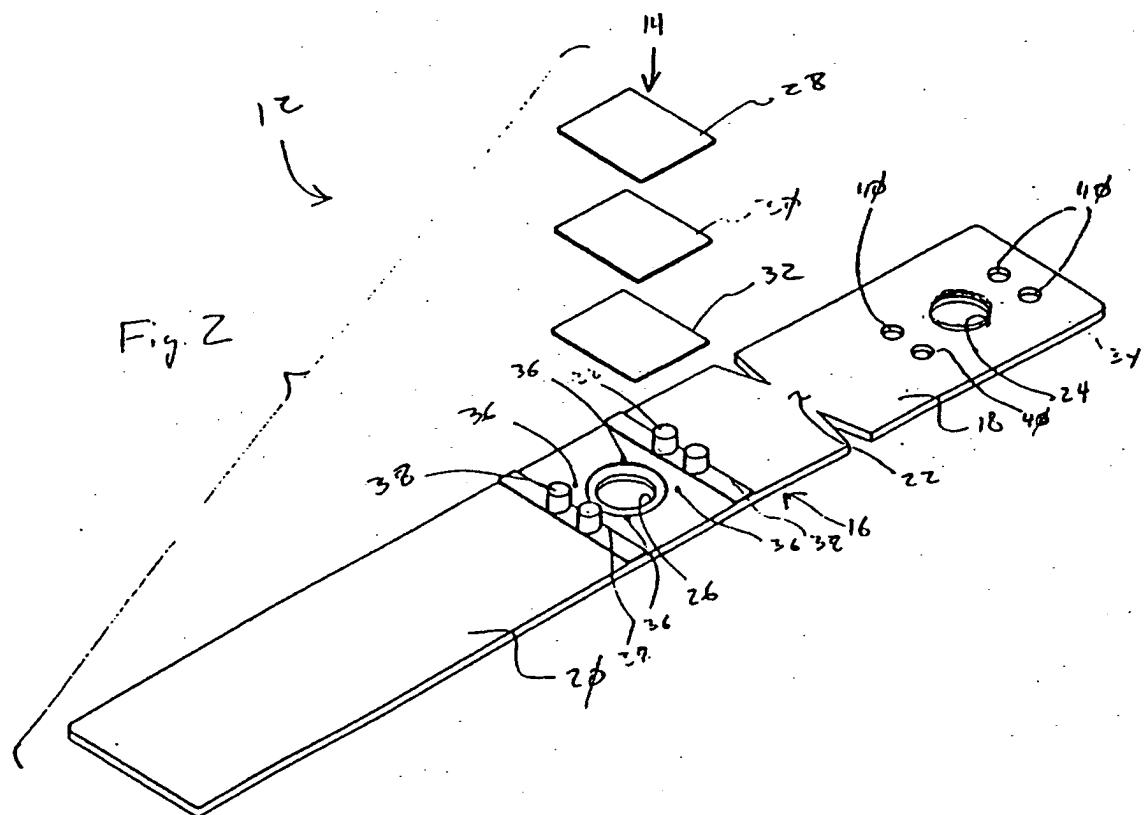


Fig. 3


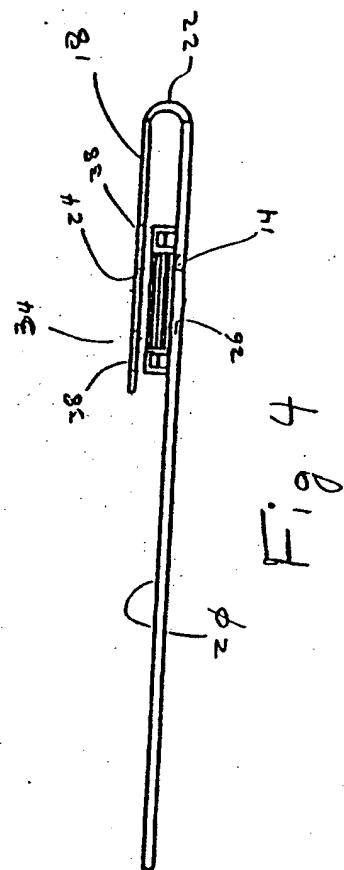


Fig. 4

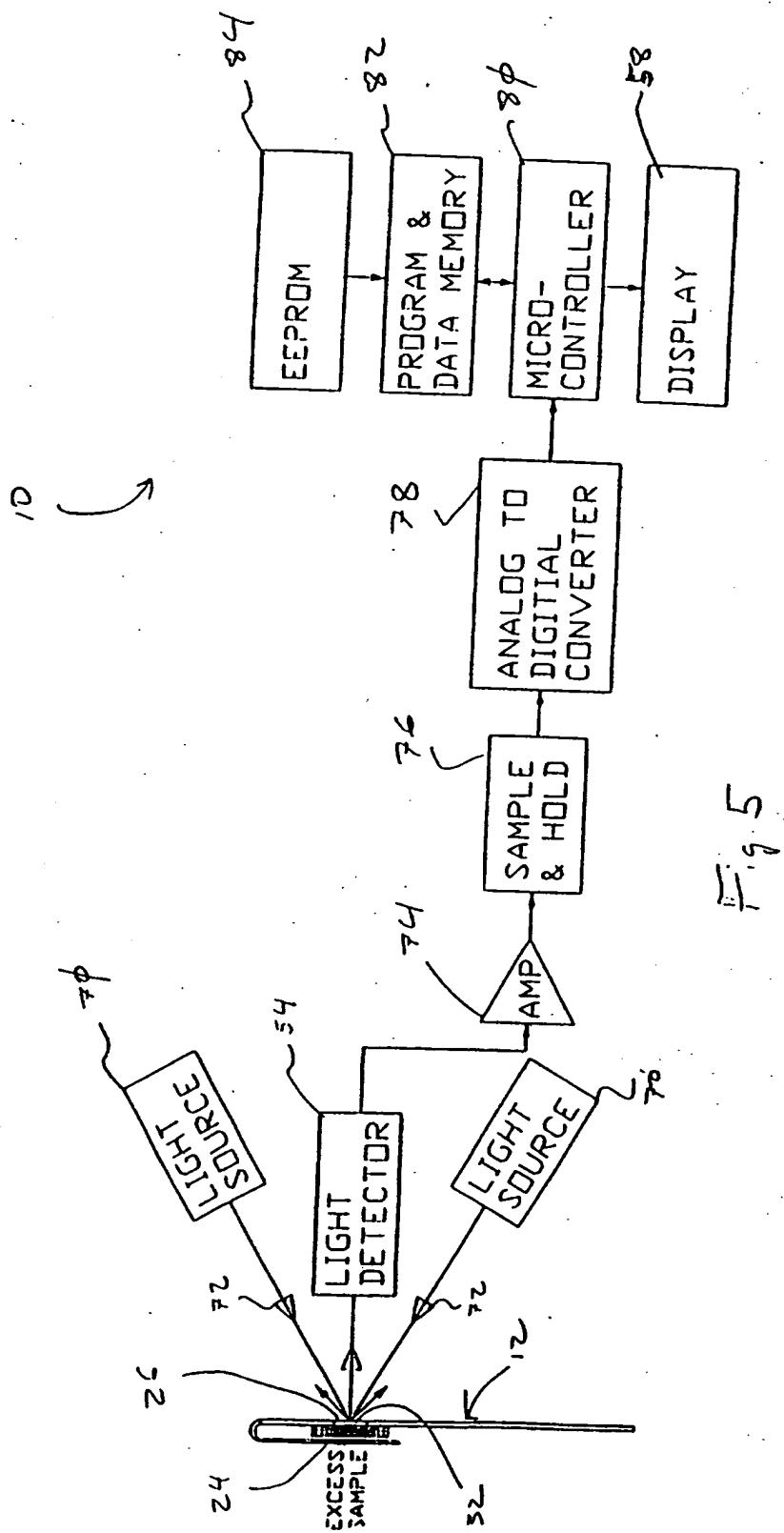

5/5

ANY REFERENCE TO  
FIGURE 6 SHALL BE  
CONSIDERED NON-EXISTENT  
(See Article 14(2))

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- BLACK BORDERS**
- IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- FADED TEXT OR DRAWING**
- BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- SKEWED/SLANTED IMAGES**
- COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- GRAY SCALE DOCUMENTS**
- LINES OR MARKS ON ORIGINAL DOCUMENT**
- REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- OTHER: \_\_\_\_\_**

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**

This Page Blank (uspto)